

# A Adam, Rouilly

AR251 - AR251-B

## Venepuncture, Cannulation, Infusion & Injection Arm

Product Leaflet



## Venepuncture, Cannulation, Infusion &amp; Injection Arm



AR251

AR251-B



This highly realistic Venepuncture and Infusion Arm Trainer is cast from life with detailed anatomy and palpable veins.

Suitable for IV injections, cannulation, blood withdrawal and infusion practice, it features silicone skin, a closed blood system, realistic flashback providing durable, practical training.

The Venepuncture and Infusion Arm Trainer features clearly defined veins in the antecubital fossa and on the back of the hand, with an exceptionally realistic finish representing a well-developed male left arm.

### Skills

- Features a flexible wrist for realism, durable silicone skin and is simple to set up, operate and maintain.
- Administering IV injections and infusions
- Intravenous cannulation on back of the hand or the forearm
- Setting up and monitoring infusions including application of dressings to secure cannula in place
- Withdrawal of blood using syringe or vacutainer

### Features



- Easy to set up
- Flexible wrist for improved realism
- Closed blood system which is clean to use
- Silicone skin moulded for improved feel, durability and realistic vein palpation
- Many procedures can be carried out before a replacement skin is required
- A realistic sensation is felt when the skin and simulated vein are penetrated
- When charged with artificial blood, a realistic 'flashback' confirms correct needle location in the vein and blood can then be withdrawn
- Continuous vein system with sealed two-bottle blood assembly reduces leakage risk
- Self-sealing veins distinguishable from blind veins by palpation, allowing repeated penetration
- The self-sealing material used for the veins, as well as the silicone skin, enhance durability
- Automatic pressure relief valve
- Replacement skins and veins are simple to fit

### Accessible Veins

- Dorsal metacarpal
- Cephalic
- Median cubital
- Basilic

### Note

21 gauge or smaller cannulae and needles must be used for practice on the arm.

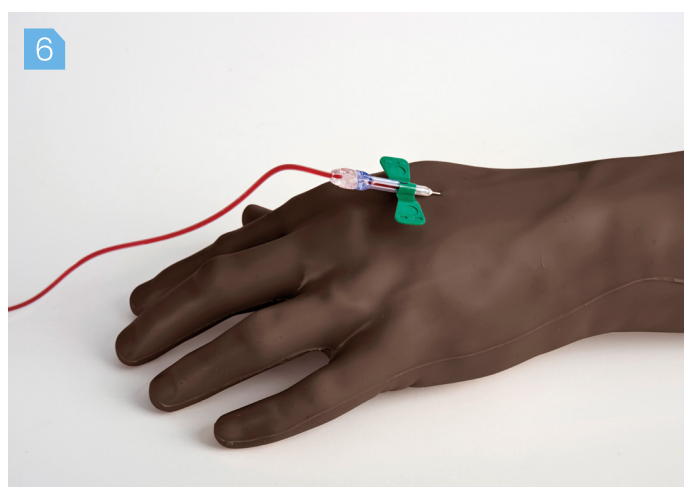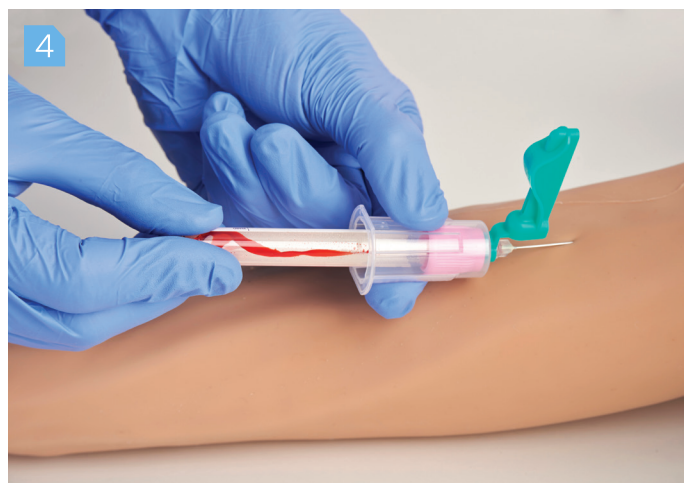

- 1 Simple Set Up
- 2 Unique Pressurised Blood System
- 3 Vein Palpation

- 4 Realistic Flashback
- 5 Cannula Insertion - Adhesive Dressing Applied
- 6 Back-of-Hand Cannula Insertion

### Includes

- Arm
- Arm Rest
- Blanking Plugs for Transportation (x2)
- Blood Bottle Assembly (x2)
- Blood Concentrate to Make Up 2.5 Litres of Artificial Blood
- Rigid Carrying Bag for Easy Storage
- Instruction Manual

### Spares & Repairs

- AR204 Blood Concentrate - Make Up 2.5 Litres
- AR253 Skin and Vein Set, Light for AR251
- AR253-B Skin and Vein Set, Dark for AR251-B
- SPAR251 Replacement Bottles, Pressure System and Arm Rest for AR251

# Safety in Training, Safety in Practice

Adam,Rouilly Limited, Castle Road, Eurolink Business Park,  
Sittingbourne, Kent, ME10 3AG, United Kingdom

## UK

+44 (0)1795 471 378  
uksales@adam-rouilly.co.uk

## Export

+44 (0)1795 411 064  
sales@adam-rouilly.co.uk

Adam,Rouilly has over 100 years of experience in providing quality medical models and simulators.

As part of our policy of continuous product development, the specifications of our products may change without prior notice.

Should you require further information, please contact our Sales Department, who will be pleased to assist you with your enquiry.

[www.adam-rouilly.co.uk](http://www.adam-rouilly.co.uk)



©Copyright 2026

**L0251 L1**